


**SHOW**  
SHORT WIRE SYSTEM

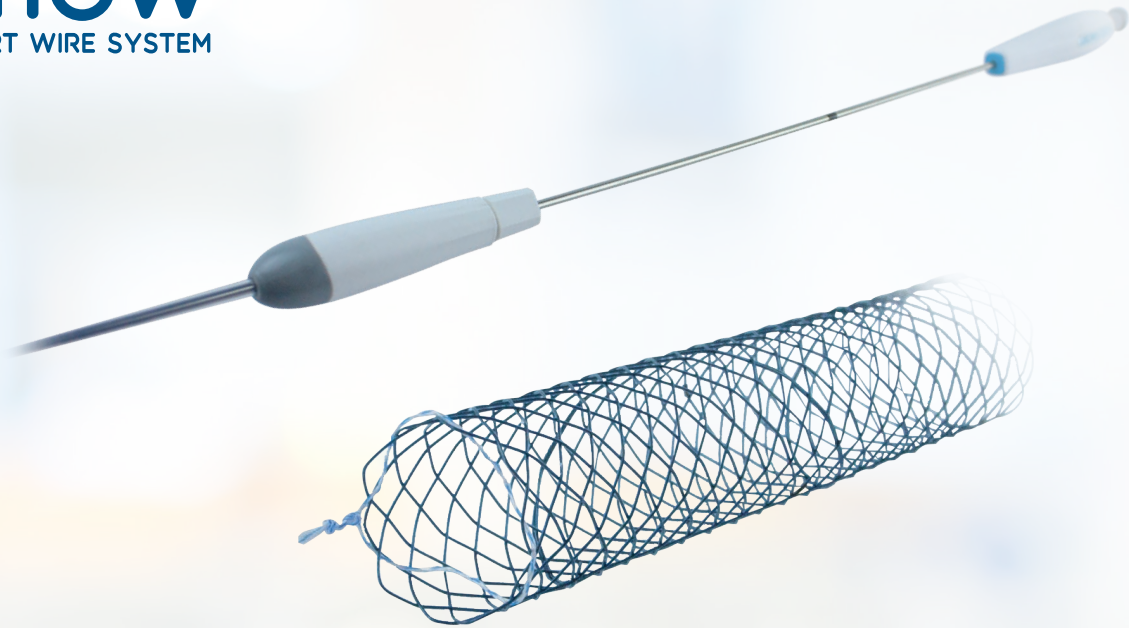

# SHOW-METAL-STENT

## INDEPENDENCE AT ITS BEST

The SHOW-Metal-Stent from MICRO-TECH is indicated for applications in the bile duct and is suitable for both short and long wire systems. Time-consuming preparation and irrigation is no longer necessary. Made of flexible nitinol wire, the stent is designed to adapt ideally to biliary conditions. The atraumatic stent ends are characterized by their soft edges. For precise placement, there are overall


10 radiopaque markers on the stent itself and three more on the introducer set. The point-of-no-return (for uncovered and partially covered stents) allows resheathing of the stent which means it can be repositioned. The guide wire exit is located at approx. 30 cm. Our range includes stents with four different lengths of 40 to 100 mm and thus offers for each operation the ideal solution.

## SPECIFIC CHARACTERISTICS

- Short and long wire compatible
- No preparation or flushing necessary
- Easy handling
- Uncovered, partially covered and fully covered
- Available in a straight form
- Automatic guide wire exit
- Point-of-no-return
- Very good visibility through tantalum markers
- Very flexible application system facilitates the placement



Extraction suture at proximal Stent-End

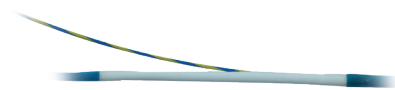

Guidewire discharge at ca. 30 cm from distal end

## SPECIFICATIONS

REF	Ø mm	Length mm	Covering mm	Ø Catheter mm	Working length mm	Repositionable	Recommended GW inch
<b>SHOW-METAL-STENTS</b>							
RST43-001-10.040	10	40	without	2.83	1900	Yes	0.035"
RST43-001-10.060	10	60	without	2.83	1900	Yes	0.035"
RST43-001-10.080	10	80	without	2.83	1900	Yes	0.035"
RST43-001-10.100	10	100	without	2.83	1900	Yes	0.035"
RST43-002-10.040	10	40	30	2.83	1900	Yes	0.035"
RST43-002-10.060	10	60	50	2.83	1900	Yes	0.035"
RST43-002-10.080	10	80	70	2.83	1900	Yes	0.035"
RST43-002-10.100	10	100	90	2.83	1900	Yes	0.035"
RST43-004-10.040	10	40	40	2.83	1900	No	0.035"
RST43-004-10.060	10	60	60	2.83	1900	No	0.035"
RST43-004-10.080	10	80	80	2.83	1900	No	0.035"
RST43-004-10.100	10	100	100	2.83	1900	No	0.035"

Package unit: 1 piece

REF	Ø inches	Total length mm	Tip length mm	Tip form	Set-up	Covering
<b>GUIDEWIRE</b>						
MTN-BM-89/26-A	0.035"	2600	65	straight	hydrophilic, radiopaque	blue-yellow
MTN-BM-89/26-A-J	0.035"	2600	65	J-form	hydrophilic, radiopaque	blue-yellow

Package unit: 2 pieces